

Module Guide for Physiotherapy Companion

Team 11, technically functional

Matthew Huynh

Cieran Diebolt

Vaisnavi Shanthamoorthy

Maham Siddiqui

Eman Ashraf

April 6, 2026

1 Revision History

Date	Version	Notes
Jan 21, 2026	1.0	Initial documentation completed.
April 6, 2026	1.1	Timeline Fixes.

2 Reference Material

This section records information for easy reference.

2.1 Abbreviations and Acronyms

symbol	description
AC	Anticipated Change
DAG	Directed Acyclic Graph
M	Module
MG	Module Guide
OS	Operating System
R	Requirement
SC	Scientific Computing
SRS	Software Requirements Specification
Physiotherapy Companion	Explanation of program name
UC	Unlikely Change

Contents

1	Revision History	i
2	Reference Material	ii
2.1	Abbreviations and Acronyms	ii
3	Introduction	1
4	Anticipated and Unlikely Changes	2
4.1	Anticipated Changes	2
4.2	Unlikely Changes	2
5	Module Hierarchy	3
6	Connection Between Requirements and Design	4
7	Module Decomposition	5
7.1	Hardware Hiding Modules	5
7.2	Behaviour-Hiding Modules	5
7.2.1	User Interface Module (M1)	6
7.2.2	Profile Activity Module (M2)	6
7.2.3	Analysis Control Flow Module (M3)	6
7.2.4	Draw Module (M4)	7
7.2.5	Account Creation Module (M5)	7
7.3	Software Decision Module	7
7.3.1	Valid User Module (M6)	7
7.3.2	User Account Module (M7)	8
7.3.3	User Activity Module (M8)	8
7.3.4	Exercise Data Module (M9)	8
7.3.5	Performance Statistics Generator Module (M10)	9
7.3.6	Metrics Module (M11)	9
7.3.7	Analyze Module (M12)	9
7.3.8	Feedback Module (M13)	10
8	Traceability Matrix	10
9	Use Hierarchy Between Modules	13
10	User Interfaces	13
11	Design of Communication Protocols	13
12	Timeline	15

List of Tables

1	Module Hierarchy	4
2	Trace Between Requirements and Modules	12
3	Trace Between Anticipated Changes and Modules	12

List of Figures

1	Use hierarchy among modules	13
2	Modules with external communications	15

3 Introduction

Decomposing a system into modules is a commonly accepted approach to developing software. A module is a work assignment for a programmer or programming team (Parnas et al., 1984). We advocate a decomposition based on the principle of information hiding (Parnas, 1972). This principle supports design for change, because the “secrets” that each module hides represent likely future changes. Design for change is valuable in SC, where modifications are frequent, especially during initial development as the solution space is explored.

Our design follows the rules laid out by Parnas et al. (1984), as follows:

- System details that are likely to change independently should be the secrets of separate modules.
- Each data structure is implemented in only one module.
- Any other program that requires information stored in a module’s data structures must obtain it by calling access programs belonging to that module.

After completing the first stage of the design, the Software Requirements Specification (SRS), the Module Guide (MG) is developed (Parnas et al., 1984). The MG specifies the modular structure of the system and is intended to allow both designers and maintainers to easily identify the parts of the software. The potential readers of this document are as follows:

- New project members: This document can be a guide for a new project member to easily understand the overall structure and quickly find the relevant modules they are searching for.
- Maintainers: The hierarchical structure of the module guide improves the maintainers’ understanding when they need to make changes to the system. It is important for a maintainer to update the relevant sections of the document after changes have been made.
- Designers: Once the module guide has been written, it can be used to check for consistency, feasibility, and flexibility. Designers can verify the system in various ways, such as consistency among modules, feasibility of the decomposition, and flexibility of the design.

The rest of the document is organized as follows. Section 4 lists the anticipated and unlikely changes of the software requirements. Section 5 summarizes the module decomposition that was constructed according to the likely changes. Section 6 specifies the connections between the software requirements and the modules. Section 7 gives a detailed description of the modules. Section 8 includes two traceability matrices. One checks the completeness of the design against the requirements provided in the SRS. The other shows the relation between anticipated changes and the modules. Section 9 describes the use relation between modules.

4 Anticipated and Unlikely Changes

This section lists possible changes to the system. According to the likeliness of the change, the possible changes are classified into two categories. Anticipated changes are listed in Section 4.1, and unlikely changes are listed in Section 4.2.

4.1 Anticipated Changes

Anticipated changes are the source of the information that is to be hidden inside the modules. Ideally, changing one of the anticipated changes will only require changing the one module that hides the associated decision. The approach adapted here is called design for change.

AC1: The user interface of the application may be changed to include additional features for accessibility. Changes include:

1. How feedback is presented to the user as the exercise is performed (i.e. written, auditory, alerts, etc.).
2. How the insights are generated to the user on their account dashboard (i.e. displayed as a line graph, columns of data, does the user have the option to select from multiple views, etc.).

AC2: The angles that are calculated to support different exercises. Currently, the application supports the 'single leg raise' exercises which is routinely used in rehabilitation for post-operative knee surgery. Future iterations of the application may include additional body parts and exercises.

AC3: Changes to the feedback/guidance given. As more requirements are met, there should be changes to encapsulate the prompts required for the motion (i.e. 'keep the torso straight as this movement is performed'). Additionally, there may be alternate prompts with varying degrees of complexity, which relates to accessibility and the patients understanding of the exercise (i.e. 'during this motion you want to engage your quadriceps as you drive your toe towards the sky 8 inches off the ground').

4.2 Unlikely Changes

The module design should be as general as possible. However, a general system is more complex. Sometimes this complexity is not necessary. Fixing some design decisions at the system architecture stage can simplify the software design. If these decision should later need to be changed, then many parts of the design will potentially need to be modified. Hence, it is not intended that these decisions will be changed.

UC1: Input/Output devices (Input: Recording/Video Upload, Output: On-Screen Feedback/Corrections).

- UC2:** The knee will remain modeled as a single-degree-of-freedom (1-DOF) hinge joint in the sagittal plane (flexion and extension only). No expansion to multi-DOF modelling or other joints is planned.
- UC3:** Users will only perform physiotherapist-approved rehabilitation exercises. The system will not automatically generate or recommend new, unverified exercises.
- UC4:** Feedback will be limited to movement-based corrections and joint angle guidance. The system will not provide medical diagnosis or replace professional physiotherapy assessment.
- UC5:** The system is not being built to integrate feedback from external sensors, wearable devices, and specialized specialized motion-capture equipment.
- UC6:** The scope will remain focused on knee rehabilitation. The system will not expand to whole-body tracking or unrelated physiotherapy movements.

5 Module Hierarchy

This section provides an overview of the module design. Modules are summarized in a hierarchy decomposed by secrets in Table 1. The modules listed below, which are leaves in the hierarchy tree, are the modules that will actually be implemented.

- M1:** User Interface Module
- M2:** Profile Activity Module
- M3:** Analysis Control Flow Module
- M4:** Draw Module
- M5:** Account Creation Module
- M6:** Account Validity Module
- M7:** User Details Module
- M8:** User Activity Module
- M9:** Exercise Data Module
- M10:** Performance Statistics Generator Module
- M11:** Metrics Module
- M12:** Analyze Module
- M13:** Feedback Module

Level 1	Level 2
Hardware-Hiding Module	
	User Interface Module
	Profile Activity Module
Behaviour-Hiding Module	Analysis Control Flow Module
	Draw Module
	Account Creation Module
	Account Validity Module
	User Details Module
	User Activity Module
Software Decision Module	Exercise Data Module
	Performance Statistics Generator Module
	Metrics Module
	Analyze Module
	Feedback Module

Table 1: Module Hierarchy

6 Connection Between Requirements and Design

The design of the system is intended to satisfy the requirements developed in the SRS. In this stage, the system is decomposed into modules. The connection between requirements and modules is listed in Table 2.

Design Decisions

Security To satisfy security and privacy requirements (FR21–FR23, SR-IR1, SR-PR1–2), the system implements end-to-end encryption for data at rest and in transit using industry-standard encryption protocols. User authentication employs session-based tokens with automatic timeout after 60 minutes of inactivity.

Real-time Performance To meet the requirement for real-time feedback (FR13–14, PR-SL1–2), the computer vision processing pipeline is optimized to provide analysis results within 1 second of detecting form deviations. This is achieved through efficient pose estimation algorithms and selective processing of critical joint angles.

Offline Capability To address reliability requirements (PR-RFT2, MS-AD1), the system implements local caching of exercise instructions and video queuing for uploads. This

ensures users can access exercise information and record sessions even without internet connectivity, with automatic synchronization when connection is restored.

Modular Exercise Library To support extensibility (PR-CR1, PR-SER1), the exercise library is designed as a modular system where new exercises and joint types can be added without modifying core analysis logic. Each exercise type has its own parameter set defining acceptable ranges and form criteria.

Two-Camera View Analysis To address the risk of 3D motion measurement with 2D cameras, the system supports analysis from the sagittal view as primary, with optional frontal view for enhanced accuracy. The analysis algorithms are designed to work with single or multiple camera angles.

7 Module Decomposition

Modules are decomposed according to the principle of “information hiding.” The *Secrets* field in a module decomposition is a brief statement of the design decision hidden by the module. The *Services* field specifies *what* the module will do without documenting *how* to do it. If the entry is *OS*, this means that the module is provided by the operating system or by standard programming language libraries. *Physiotherapy Companion* means the module will be implemented by the Physiotherapy Companion software. Only the leaf modules in the hierarchy have to be implemented.

7.1 Hardware Hiding Modules

Secrets: The data structure and algorithm used to implement the virtual hardware.

Services: Serves as a virtual hardware used by the rest of the system. This module provides the interface between the hardware and the software, including camera access, screen display, audio output, and device sensors. The system can use it to capture video, display outputs, or accept touch inputs.

Implemented By: Android OS

7.2 Behaviour-Hiding Modules

Secrets: The contents of the required behaviours.

Services: Includes programs that provide externally visible behaviour of the system as specified in the software requirements specification (SRS) documents. This module serves as a communication layer between the hardware-hiding module and the software decision module. The programs in this module will need to change if there are changes in the SRS.

Implemented By: Physiotherapy Companion

7.2.1 User Interface Module (M1)

Secrets: The user interface layout and interaction logic for patient and physiotherapist dashboards.

Services: Displays patient exercise history, progress charts, and upcoming exercises. Shows PT patient list with status indicators and alerts. Provides access to video review and report viewing. Enables PT comments on patient sessions. Presents analytics and visualizations of exercise data.

Implemented By: Physiotherapy Companion

Type of Module: Abstract Object

Related Requirements: FR18, FR19, FR24, FR25

7.2.2 Profile Activity Module (M2)

Secrets: The data structures and business logic for managing patient-specific exercise programs, including customization and progress tracking.

Services: Allows physiotherapists to create and customize exercise programs for individual patients, including sets, reps, rest times, and weekly limits. Tracks patient exercise completion and progress over time. Enforces PT-defined limitations on exercise frequency and intensity.

Implemented By: Physiotherapy Companion

Type of Module: Abstract Object

Related Requirements: FR6, FR9, FR10, FR11, FR22, FR27, SR-AR3

7.2.3 Analysis Control Flow Module (M3)

Secrets: The implementation of video recording, including camera configuration, frame rate management, and recording controls.

Services: Captures video of users performing exercises at minimum 720p resolution and 30 fps. Performs pre-capture environment checks (lighting, framing). Provides recording controls (start, stop, retry). Ensures video meets minimum quality standards before proceeding to analysis.

Implemented By: Physiotherapy Companion

Type of Module: Abstract Object

Related Requirements: FR7, FR9, FR13, FR15, FR26, UHR-EOU3

7.2.4 Draw Module (M4)

Secrets: The graphical representation of landmarks and objects of interest within the video provided by the user.

Services: Graphically shows the user the system's understanding of their movements by highlighting landmarks such as their ankle, knee, and hip. Additionally, draws lines between these landmarks to communicate how the system interprets limb movement.

Implemented By: Physiotherapy Companion

Type of Module: Library

Related Requirements: FR5, FR14, LFR-ST1, UHR-EOU2, UHR-EOU5, UHR-EOU6

7.2.5 Account Creation Module (M5)

Secrets: The logic for creating and storing both PTs and patient accounts, and creating physiotherapist invite codes.

Services: Handles new user account creations for both patients and PTs with input information. Also generates unique invite codes for PT accounts upon creation for linking patient accounts. Stores user account information securely in a remote database.

Implemented By: Physiotherapy Companion

Type of Module: Library Module

Related Requirements: FR1, FR2

7.3 Software Decision Module

Secrets: The design decision based on mathematical theorems, physical facts, or programming considerations. The secrets of this module are not described in the SRS.

Services: Includes data structure and algorithms used in the system that do not provide direct interaction with the user.

Implemented By: Physiotherapy Companion

7.3.1 Valid User Module (M6)

Secrets: The validation and authorization logic for all user credentials (patient and PT), invite codes, and physiotherapist identification numbers.

Services: Checks input user information to ensure unique entries into the remote database. Ensures valid invite code during patient account creation to link new account to relevant PT account. Ensures physiotherapist identification numbers are valid through the publicly accessible records during PT account creation.

Implemented By: Physiotherapy Companion

Type of Module: Library Module

Related Requirements: FR1, FR3

7.3.2 User Account Module (M7)

Secrets: The authorization logic for existing user actions, and the retrieval and modification logic for user information in the remote database.

Services: User logins (both patient and PT) authorized by checking the remote database for matching email and password. Enables physiotherapists to delete linked patient accounts for housekeeping and service limitation purposes. Additionally, retrieves all relevant user account information from the remote database for other module usage.

Implemented By: Physiotherapy Companion

Type of Module: Data Access Module

Related Requirements: FR1, FR2, FR4, FR18, FR21, FR22, FR23, SR-AR1, SR-AR3, SR-IR1, SR-PR1, SR-PR2

7.3.3 User Activity Module (M8)

Secrets: The structure and storage logic for users historic activity data.

Services: Formats, links, and stores user activity data in the remote database in a table for historic activity summaries. Formats input from the Analyze Module.

Implemented By: Physiotherapy Companion

Type of Module: Abstract Object

Related Requirements: FR21, FR25, SR-IR1

7.3.4 Exercise Data Module (M9)

Secrets: The structure and storage logic for user relevant exercise data.

Services: Provides exercise information to users on demand. Maintains a library of knee recovery relevant exercises along with exercise instructions. Allows user specific modification of prescribed exercises by linked PT account.

Implemented By: Physiotherapy Companion

Type of Module: Abstract Object

Related Requirements: FR5, FR6, FR13, FR21, MS-AD1, SR-IR1, SR-PR1

7.3.5 Performance Statistics Generator Module (M10)

Secrets: The format and generation logic for exercise session reports and progress summaries.

Services: Generates comprehensive session reports including completion rate, form score, pain level, and patient feedback. Creates progress summaries showing improvement over time. Generates alerts for repeated incorrect attempts or concerning patterns. Formats reports for PT review.

Implemented By: Physiotherapy Companion

Type of Module: Library

Related Requirements: FR17, FR19, FR20, FR25

7.3.6 Metrics Module (M11)

Secrets: The specific computer vision algorithms and implementation details for pose detection, joint tracking, and body landmark identification.

Services: Processes video frames to detect and track body landmarks using pose estimation. Identifies key joints relevant to knee exercises (hip, knee, ankle, foot). Extracts coordinate data for joint positions across video frames. Handles varying lighting conditions and camera angles. Implements the core CV library integration (OpenCV, MediaPipe).

Implemented By: Physiotherapy Companion

Type of Module: Abstract Data Type

Related Requirements: FR8, FR12

7.3.7 Analyze Module (M12)

Secrets: The mathematical models and algorithms for calculating joint angles, range of motion, and determining form correctness based on exercise parameters. Also the mechanism for storing video files through local caching.

Services: Calculates joint angles from landmark coordinates with $\pm 5^\circ$ accuracy. Measures range of motion with $\pm 10\%$ accuracy. Compares measured values against exercise-specific acceptable ranges. Detects form deviations and classifies severity. Identifies dangerous movements requiring immediate halt. Measures time-under-tension and exercise duration. Also stores recorded exercise videos securely (encrypted at rest per SR-IR1).

Implemented By: Physiotherapy Companion

Type of Module: Abstract Data Type

Related Requirements: FR12, FR13, FR14, FR21, SR-IR1, PR-RFT2, PR-PAR1, PR-PAR2, PR-SCR1, MS-MNT1

7.3.8 Feedback Module (M13)

Secrets: The rules and logic for converting analysis results into user-friendly feedback messages and instructions.

Services: Generates real-time corrective feedback messages based on form analysis results. Provides step-by-step exercise guidance. Creates descriptive messages for form deviations specifying required adjustments. Triggers safety warnings for dangerous movements.

Implemented By: Physiotherapy Companion

Type of Module: Library

Related Requirements: FR13, FR14, FR16, PR-SCR1

8 Traceability Matrix

This section shows two traceability matrices: between the modules and the requirements and between the modules and the anticipated changes.

Req.	Modules
FR1	M5, M6, M7
FR2	M5, M7
FR3	M6
FR4	M7
FR5	M4, M9
FR6	M2, M9, M10

FR7	M3
FR8	M3, M11
FR9	M2, M3
FR10	M2
FR11	M2
FR12	M11, M12
FR13	M3, M9, M12, M13
FR14	M4, M12, M13
FR15	M3
FR16	M13
FR17	M10
FR18	M1, M7
FR19	M1, M10
FR20	M10
FR21	M7, M8, M9, M12
FR22	M2, M7
FR23	M7
FR24	M1
FR25	M1, M8, M10
FR26	M3
FR27	M2
LFR-ST1	M4
MS-AD1	M9
MS-MNT1	M12
PR-PAR1	M12
PR-PAR2	M12
PR-RFT2	M12
PR-SCR1	M12, M13
SR-AR1	M7
SR-AR2	M6
SR-AR3	M2, M7
SR-IR1	M7, M8, M9, M12
SR-PR1	M7, M9
SR-PR2	M7
UHR-EOU2	M4

UHR-EOU3	M3
UHR-EOU5	M4
UHR-EOU6	M4

Table 2: Trace Between Requirements and Modules

AC	Modules
AC1	M2, M3, M10, M13
AC2	M2, M3, M8, M9, M11, M12
AC3	M3, M11, M12, M13

Table 3: Trace Between Anticipated Changes and Modules

9 Use Hierarchy Between Modules

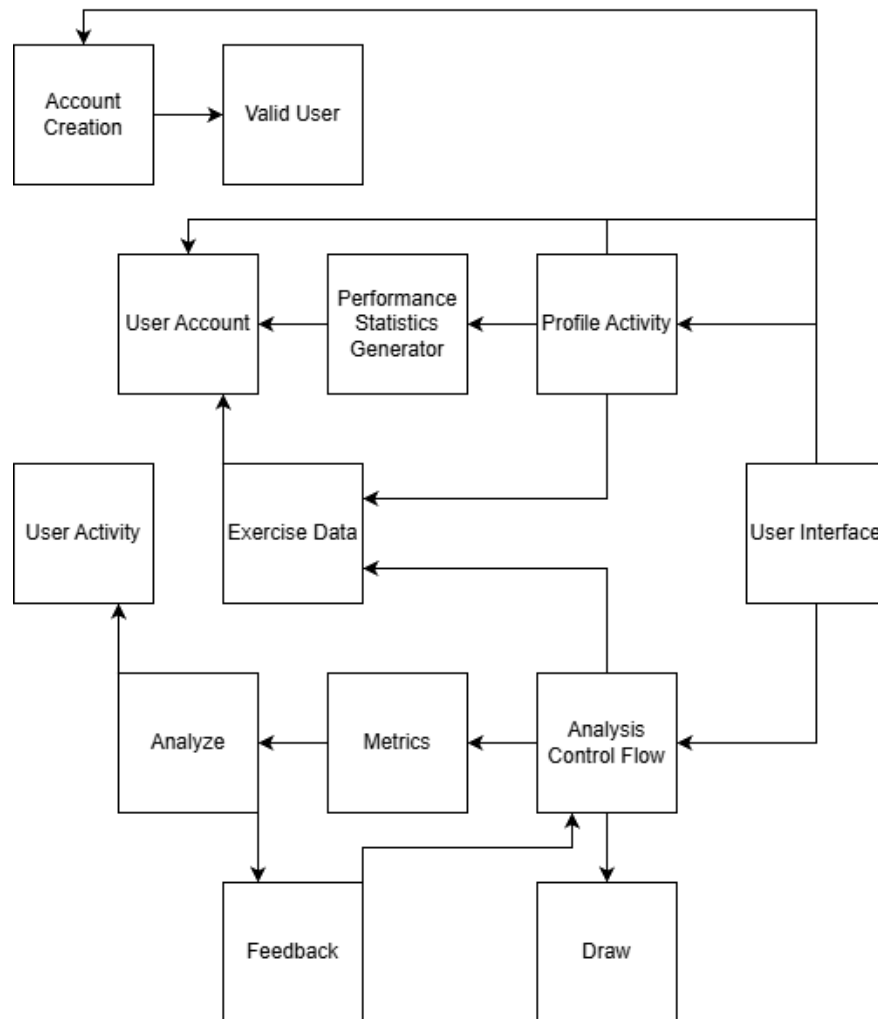


Figure 1: Use hierarchy among modules

10 User Interfaces

The proposed user interface for the system will attempt to model [this Figma prototype](#).

11 Design of Communication Protocols

Hardware intake is limited to the video feed. This feed, provided by the respective system’s hardware layer, is requested and used by the Analysis Control Flow Module. This module also overlays respective data returned to it from the Metrics Module. Additionally,

the metrics module makes API calls to OpenCV and MediaPipe API to determine necessary information such as joint angle, movement speed, exercise duration, and number of repetitions.

The Metrics Module also passes this relevant information to the Analyze Module in real-time where the numerical returns are inspected to ensure exercises are being performed within bounds. The Analyze Module makes necessary calls to the Feedback Module if the user should be told to change form, speed, or other such exercise relevant information. The Analyze Module also passes generalized exercise information to the User Activity Module for storing exercise overviews tied to user accounts.

The Account Creation, Valid User, User Account, User Activity, and Exercise Data Modules interact with the external database hosted by Firebase.

The Account Creation Module uploads new account information to a table kept for user information which includes details relevant for login and account validation.

The Valid User Module checks the user information table to ensure an account doesn't exist with identical information and that the entered invitation code is valid to an existing physiotherapist account.

The User Account Module both uploads and retrieves information from the user information table for usage on the main account page and for linking relevant user activities and exercise data.

The User Activity Module uploads information to an activity tracking table in the database. This table is used to show both the user and the linked physiotherapist the results and overall progress on prescribed exercises.

The Exercise Data Module both uploads and retrieves information from a table for specific exercise alterations. This information is editable by physiotherapist accounts and is referenced by the Analysis Control Flow Module during exercise recording.

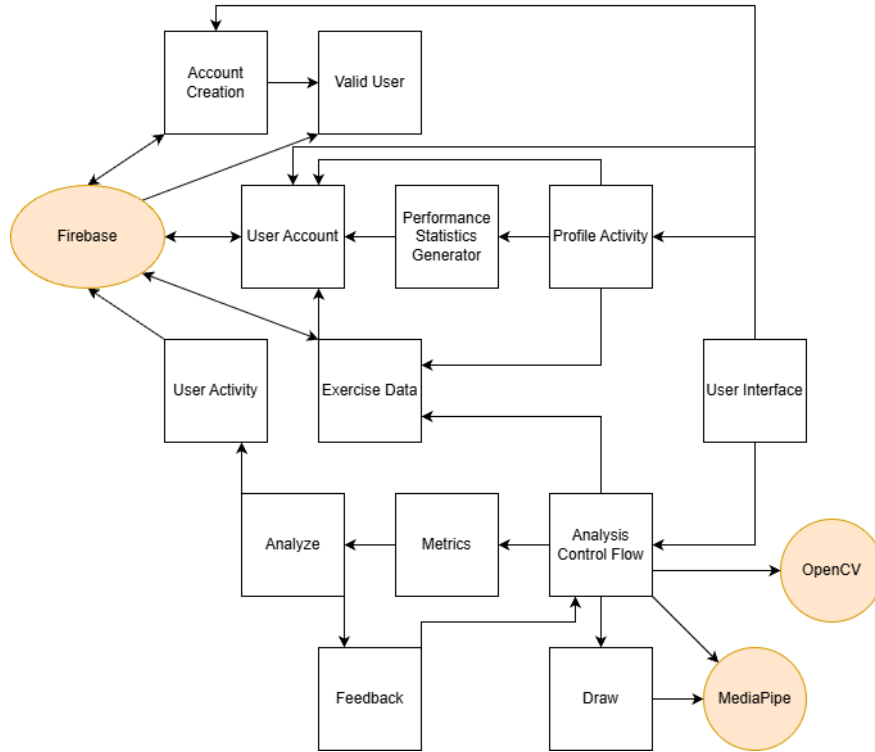


Figure 2: Modules with external communications

12 Timeline

For the remainder of the term, the project timeline will follow the official course deliverables and our internal task assignments. A detailed task breakdown is maintained on our [GitHub Project Board](#). The high-level schedule is listed below:

- **November 14th, 2025: Analysis Control Flow Module** (Responsible: Matthew, Vaisnavi)
 - Video recording capabilities fully implemented and ready for POC.
- **November 14th, 2025: Draw Module** (Responsible: Matthew)
 - Landmark placement and line drawing capabilities fully implemented and ready for POC.
- **November 14th, 2025: Metrics Module** (Responsible: Matthew, Vaisnavi)
 - Numerical data extraction from video feed fully implemented and ready for POC.
- **November 14th, 2025: Analyze Module** (Responsible: Vaisnavi)

- Exercise analysis, feedback generation, and data parsing capabilities fully implemented and ready for POC.
- **November 16th, 2025: POC Testing** (Responsible: Entire Team)
 - Run throughs and debugging of the POC to ensure everything is working as expected before the POC demonstration.
- **During November, 2025: Design Doc Rev0 and further development** (Responsible: Eman, Cieran, Maham)
 - Work on the remaining modules (see dates below), and refine design documents from -1 to 0 for deliverable in January.
- **January 23rd, 2026: Design Doc Rev0** (Responsible: Entire Team)
 - Complete MG and MIS documents with all modules fully described and connected to requirements, and with the use hierarchy fully defined.
- **January 26th, 2026: Account Creation Module** (Responsible: Maham)
 - Complete implementation of account creation capabilities, including user input handling and database communication.
- **January 26th, 2026: UI Module** (Responsible: Vaisnavi)
 - React frontend development of the user interface with basic navigational capabilities, and respective places for linking backend data to be displayed.
- **January 30th, 2026: Valid User Module** (Responsible: Matthew, Vaisnavi)
 - Complete camera and pose detection reliability with MediaPipe.

December to January

- **Work on project as a whole in preposition for Revision 0 Demo** (Responsible: Entire Team)
 - Work on entire project, start building and revising the remaining modules, to begin putting together all the pieces

February

- **Rev0 Presentations** (Responsible: Entire Team)
 - Implement feedback from the professor and TA in the previous deliverables.
 - Refactor code to ensure it is functional, stable, and ready for the presentation.
 - Finalize demonstration plan and speaking roles.

Late March – Early April

- **Final Demonstration (Revision 1)** (Responsible: Entire Team)
 - Implement improvements from TA feedback.
 - Enhance UI transparency (tracking visualization, angles, notifications).
 - Prepare live demo and video recording.
- **EXPO Demonstration** (Responsible: Entire Team)
 - Build poster and prepare speaking roles.
 - Conduct demo rehearsals and ensure hardware stability.

April – Final Week

- **Final Documentation (Revision 1)** (Responsible: Entire Team)
 - Submit Problem Statement, Development Plan, POC Plan, Requirements, Hazard Analysis, Design Document, V&V Plan, V&V Report, Extra Documentation, and Source Code.
 - Include GitHub repository link, README, demo video, and traceability matrices.

References

- David L. Parnas. On the criteria to be used in decomposing systems into modules. *Comm. ACM*, 15(2):1053–1058, December 1972.
- D.L. Parnas, P.C. Clement, and D. M. Weiss. The modular structure of complex systems. In *International Conference on Software Engineering*, pages 408–419, 1984.